



Product Catalogue

Quality Supplies for Hearing
Healthcare Professionals

Sanibel
Supply

Valid from 1 August 2024

Contents

- 05** Disposable ear tips
- 09** Probe tips
- 11** Supplies for insert earphones
- 13** ABR and Newborn
- 15** Disposable electrodes
- 19** Electrode cables
- 20** Patient preparation
- 22** Infection control
- 24** REM probe tubes for HA fitting
- 25** Hearing aid fitting and supplies
- 32** Clinical tools
- 34** Hearing aid maintenance and hearing care
- 36** Cerumen Management accessories
- 42** Balance supplies
- 43** Paper
- 45** Otoscopy supplies
- 47** Accessories



Why use ear tips for single use only?

Sanibel Supply® recommends using disposable ear tips for single use only to stop ear infections from spreading

At Sanibel Supply® we think that most audiologists would agree with us that it is crucial to avoid inadvertent transmission of disease within the audiology clinic. The good news is that it is quite easy to ensure that ear infections are not spread from one patient to another; use a new ear tip for each patient, and you will significantly reduce or eliminate the incidences of cross contamination.

Cerumen may be an infectious agent

Cerumen is not considered an infectious agent unless contaminated with blood or mucus, and due to its color and viscosity, visual detection of blood or ear discharge contaminants may be difficult. Therefore, you should treat cerumen as a potentially infectious agent. It is especially important to be mindful of using new tips when testing children and those with weaker immune systems, such as the frail and elderly as they are more susceptible to infections.

Are disposable ear tips expensive?

Sanibel Supply® offers volume discounts on ear tips. This will often offset the time and cost involved in the ear tip cleaning process.

Be more cost effective in your clinic

By using a new ear tip pr. every patient, you reach a seal much more reliable than with re-used ear tips. How often have you not been fumbling for minutes before you get a seal with a patient? Using the right ear tip and a fresh one pr. patient ensures you a reliable fit with every patient, meaning you can calculate less minutes for testing a patient and in most cases allow you to have an extra patient pr. day without having stressed out personnel. The cost of using one ear tip pr. patient is far less than what you earn on your next patient, plus you can promote your clinic as taking patient safety and hygiene seriously.

6 strong reasons why you should use a new ear tip after each patient:

- 1**
You reduce the risk of cross contamination between your patients
Unfortunately, it is not possible to be 100% certain that complete decontamination has taken place when cleaning ear tips.
- 2**
You minimize the risk of non-sealing ear tips
Ear tips degrade over time, which can make it more difficult to attain a good seal, hence leading to less accurate results.
- 3**
Your patients avoid allergic reactions from cleaning agent residue
Residue from cleaning agents may cause allergic reactions in some patients.
- 4**
Your patients will appreciate it
Patients feel safer and are more comfortable knowing that new ear tips are being used.
- 5**
It's so much easier
Cleaning ear tips takes time, so why not give your staff the time to do what really matters?
- 6**
It's so much faster
A fresh ear tip creates a seal quicker and more reliable, so you do not waste time performing a standard procedure.

List of abbreviations

ABR	Auditory Brainstem Response	HIF	Hearing Instrument Fitting
AUD	Audiometry	IMP	Impedance
BAL	Balance	OAE	Otoacoustic Emissions



Single patient use only!



Product is not made with natural rubber latex




























Selecting the correct ear tip

Selecting the correct size of ear tip takes practice

In general, flanged or mushroom-shaped ear tips will be used for most testing (Impedance, DPOAE, TEOAE and Automated ABR).

Sometimes, an umbrella-shaped ear tip will be used for performing a quick impedance screening.

The table can be used as a guide to help you select the appropriate ear tip size.

Age	Ear tips inserted against ear canal	Ear tips inserted against entrance of ear canal			Ear tips for Echo-Screen®	Ear tips for AccuScreen®
0.5 y/o			Green, pediatric		 4.0 mm, ANA	 3.7 mm, AOT
					 4.5 mm, ANA	 4.0 mm, AOT
			3-5 mm, flanged	 3-6 mm, green	 6 mm, green	 4.5 mm, AOT
	 7.5 mm, umbrella	 4-7 mm, flanged	 6 mm, yellow		 infant tree, ANA	 5.0 mm, AOT
	 9 mm, umbrella	 5-8 mm, flanged			 mushroom, ANA	 4-7 mm, AOT
10 y/o	 10 mm, umbrella	 7 mm, mushroom	 7 mm, red	 7 mm, yellow		
	 11 mm, umbrella	 8 mm, mushroom	 8 mm, blue	 8 mm, red		
	 13 mm, umbrella	 9 mm, mushroom	 9 mm, green	 9 mm, blue		
		 10 mm, mushroom	 10 mm, yellow	 10 mm, green		
	 15 mm, umbrella	 11 mm, mushroom	 11 mm, red	 11 mm, yellow		
Adults		 12 mm, mushroom	 12 mm, blue	 12 mm, red		
	 19 mm, umbrella	 13 mm, mushroom	 13 mm, green	 13 mm, blue		
		 14 mm, mushroom	 14 mm, yellow	 14 mm, green		
		 15 mm, mushroom	 16 mm, red	 16 mm, yellow		
	 22 mm, umbrella	 19 mm, mushroom	 18.5 mm, blue			
						 Medical grade silicone rubber Biocompatible

Disposable ear tips

Disposable Sanibel™ ADI ear tips

For use with the following acoustic immittance instruments and otoacoustic emission systems with CLEAR probe tip: **Amplivox**: all Otowave versions. **Interacoustics**: Titan, AT235, AA220, Eclipse, MT10, Sera, OtoRead II. **MAICO**: easyTymp, touchTymp MI 22, MI 23, MI 24, MI 26, MI 34, MI 44, easyScreen, EroScan Plus.

Alternative to AKR (KR) / IA series: Old MAICO MI-series + RaceCar black probe versions together with adaptor kit (8500430). **GSI**: Novus, TympStar Pro, Corti and Allegro. **Natus**: Zodiac (not 901) and AccuScreen.

\$105



ADI Assortment Box

Part no.	Description
8511445	Various sizes for tympanometry and OAE

\$70



Flanged infant ear tips

Part no.	Description	Unit
8502838	Pediatric (green)	Bag/100 pcs
8012964	3-5 mm (red)	Bag/100 pcs
8012966	4-7 mm (blue)	Bag/100 pcs
8012968	5-8 mm (yellow)	Bag/100 pcs

2

LATEX

ABR

IMP

OAE

\$70



Mushroom ear tips

Part no.	Description	Unit
8013002	7 mm (blue)	Bag/100 pcs
8013004	8 mm (yellow)	Bag/100 pcs
8012970	9 mm (green)	Bag/100 pcs
8012972	10 mm (red)	Bag/100 pcs
8012974	11 mm (blue)	Bag/100 pcs
8012976	12 mm (yellow)	Bag/100 pcs
8012978	13 mm (green)	Bag/100 pcs
8012980	14 mm (red)	Bag/100 pcs
8012982	15 mm (blue)	Bag/100 pcs
8012984	19 mm (yellow)	Bag/100 pcs

2

LATEX

ABR

IMP

OAE

\$70



Umbrella ear tips

Part no.	Description	Unit
8012986	7.5 mm (green)	Bag/100 pcs
8012988	9 mm (red)	Bag/100 pcs
8012990	10 mm (blue)	Bag/100 pcs
8012992	11 mm (yellow)	Bag/100 pcs
8012994	13 mm (green)	Bag/100 pcs
8012996	15 mm (red)	Bag/100 pcs
8012998	19 mm (blue)	Bag/100 pcs
8013000	22 mm (yellow)	Bag/100 pcs

2

LATEX

IMP

\$70

Single use ear tips for immittance instruments & OAE screeners

For use with many acoustic immittance instruments and otoacoustic emission systems, including the **MAICO** EroScan, **Interacoustics** OtoRead and **Grason-Stadler** TympStar.



2

LATEX

ABR

OAE

Infant ear tips

Part no.	Description	Unit
8013006	3 mm (yellow)	Bag/100 pcs
8002035	4 mm flanged (red)	Bag/100 pcs

Disposable ear tips

IA series ear tips

For use with many acoustic immittance instruments and otoacoustic emission systems, including: **GSI**: TympStar Pro, Allegro, Novus and Corti. **Interacoustics** (clear probe): AT235, AA222 and Titan and Otoread II, MT10. **MAICO**: touchTympt RaceCar, touchTympt MI24, MI26, MI34, MI36, easyTympt. **Madsen/Otometrics**: Zodiac (new), Accuscreen. **Path Medical**: all Sentiero products.

\$90



Mushroom ear tips

Part no.	Description	Unit
8520348	3-6 mm (green)	Bag/100 pcs
8520349	6 mm (yellow)	Bag/100 pcs
8520350	7 mm (red)	Bag/100 pcs
8520351	8 mm (blue)	Bag/100 pcs
8520352	9 mm (green)	Bag/100 pcs
8520353	10 mm (yellow)	Bag/100 pcs
8520354	11 mm (red)	Bag/100 pcs
8520355	12 mm (blue)	Bag/100 pcs
8520356	13 mm (green)	Bag/100 pcs
8520357	14 mm (yellow)	Bag/100 pcs
8520358	16 mm (red)	Bag/100 pcs
8520359	18.5 mm (blue)	Bag/100 pcs

2

LATEX

ABR

IMP

OAE

Disposable Sanibel™ AOT and ANA ear tips

\$48



ANA ear tips for Echo-Screen®

For use with **Natus**: Echo-Screen®.
Madsen/Otometrics: AccuScreen® Classic.
Interacoustics: Eclipse with pink probe tip

Part no.	Description	Unit
8106363	4.0 mm (blue)	Bag/50 pcs
8106362	4.5 mm (red)	Bag/50 pcs
8106364	5.0 mm (green)	Bag/50 pcs
8106366	Infant tree (yellow)	Bag/50 pcs
8106365	Mushroom (green)	Bag/50 pcs

2

LATEX

ABR

OAE

\$48



AOT ear tips

For use with **Madsen/Otometrics** AccuScreen®, Biologic Bio-logic® AuDX®. **Path Medical**: all Sentiero products

Part no.	Description	Unit
8106367	3.7 mm (yellow)	Bag/50 pcs
8106368	4.0 mm (light green)	Bag/50 pcs
8106369	4.5 mm (light pink)	Bag/50 pcs
8106370	5.0 mm (light blue)	Bag/50 pcs
8106371	4-7 mm (transparent)	Bag/50 pcs

2

LATEX

ABR

OAE

AKR mushroom ear tips

Use with the following acoustic immittance instruments and otoacoustic emission systems with BLACK probe tip from **Interacoustics**: AT235, AA220, Eclipse, MT10, Otoread (I), **MAICO**: Old RaceCar, old MI22-36, MI 22, MI 23, MI 24, MI 26, MI 34, MI 44. **Amplivox**: Otowave (all with white probe). **Madsen/Otometrics**: Capella 2, Zodiac 901, Biologic AuDX Pro (+Flex). **Neurosoft**: Neuro-Audio, **Oscilla**: TSM 300/400 + T800 series. All **GSI** tympanometers from GSI 20 - 39 + Tympstar, Allegro.

\$81



Mushroom ear tips

Part no.	Description	Unit
8500570	3-6 mm (yellow)	Bag/100 pcs
8107510	4 mm (red)	Bag/100 pcs
8107511	5 mm (blue)	Bag/100 pcs
8500560	6 mm (green)	Bag/100 pcs
8500550	7 mm (yellow)	Bag/100 pcs
8500540	8 mm (red)	Bag/100 pcs
8500530	9 mm (blue)	Bag/100 pcs
8500520	10 mm (green)	Bag/100 pcs
8500510	11 mm (yellow)	Bag/100 pcs
8500500	12 mm (red)	Bag/100 pcs
8500490	13 mm (blue)	Bag/100 pcs
8107520	14 mm (green)	Bag/100 pcs
8500480	16 mm (yellow)	Bag/100 pcs

2

LATEX

ABR

OAE

\$277



Eartip storage cabinet

12 small drawers, 3 medium drawers, 1 large drawer.
 Can be hung on the wall. H 28.2 cm x W 30.6 cm x L 15 cm

Part no.	Unit
8100665	1

IT Select Series Ear Tips

For use with Inventis: Flute, Clarinet, Viola. Resonance: r16m, r26m, r36m, r15c, r25c, r35c

\$80



IT Select

Part no.	Description	Unit
8522402	8mm (green)	Bag/100 pcs
8522403	10mm (red)	Bag/100 pcs
8522404	12mm (blue)	Bag/100 pcs
8522405	14mm (green)	Bag/100 pcs
8522406	16mm (red)	Bag/100 pcs

Disposable ear tips

Grason & Associates OT & MO select series ear tips

\$81



OT select ear tips w/ triangle holes

Ear tips with triangle holes for use with
Otodynamics UGD, SND, SGD probes

Part no.	Description	Unit
8102644	4.8 mm (blue)	Bag/100 pcs
8102645	6 mm (green)	Bag/100 pcs
8102646	8 mm (red)	Bag/100 pcs
8102647	10 mm (blue)	Bag/100 pcs
8102648	12 mm (green)	Bag/100 pcs
8102649	14 mm (red)	Bag/100 pcs
8102650	16 mm (blue)	Bag/100 pcs

2

LATEX

ABR

OAE



OT select ear tips w/two probe holes

Ear tips with two probe holes for use with
Otodynamics ECP, UGS, SNS, SGS probes.

Part no.	Description	Unit
8102632	3.5 mm (blue)	Bag/100 pcs
8102633	4.5 mm (blue)	Bag/100 pcs
8102634	5.5 mm (blue)	Bag/100 pcs
8102635	6.5 mm (blue)	Bag/100 pcs
8102636	4.8 mm (blue)	Bag/100 pcs
8102637	6 mm (green)	Bag/100 pcs
8102638	8 mm (red)	Bag/100 pcs
8102639	10 mm (blue)	Bag/100 pcs
8102640	12 mm (green)	Bag/100 pcs
8102641	14 mm (red)	Bag/100 pcs
8102642	16 mm (blue)	Bag/100 pcs

\$81

2

LATEX

ABR

OAE



MO select ear tips for Madsen OTOflex

For **MADSEN**® OTOflex™.

Part no.	Description	Unit
8102653	3-6 mm (yellow)	Bag/100 pcs
8102654	6 mm (red)	Bag/100 pcs
8102655	7 mm (blue)	Bag/100 pcs
8102656	8 mm (green)	Bag/100 pcs
8102657	9 mm (yellow)	Bag/100 pcs
8102658	10 mm (red)	Bag/100 pcs
8102659	11 mm (blue)	Bag/100 pcs
8102660	12 mm (green)	Bag/100 pcs
8102661	13 mm (yellow)	Bag/100 pcs
8102662	14 mm (red)	Bag/100 pcs
8102663	16 mm (blue)	Bag/100 pcs
8102664	18 mm (green)	Bag/100 pcs

\$80

2

LATEX

ABR

OAE

Ear tips for OAE screeners

For use with **Biologic AuDx**, **Welch Allyn** OAE screener
and **Natus** Echo-Screen® III.

\$220



OAE ear tips, AuDx


Part no.	Description	Unit
8107885	AuDx Preemie 3 mm	Bag/100 pcs
8107889	AuDx Tree 3.5-5mm	Bag/100 pcs
8107887	AuDx Infant 6 mm	Bag/100 pcs
8102681	Adult OAE foam	Bag/100 pcs
8102682	Pediatric OAE foam	Bag/100 pcs
8103123	Jumbo OAE foam	Bag/100 pcs

2

LATEX

ABR

OAE

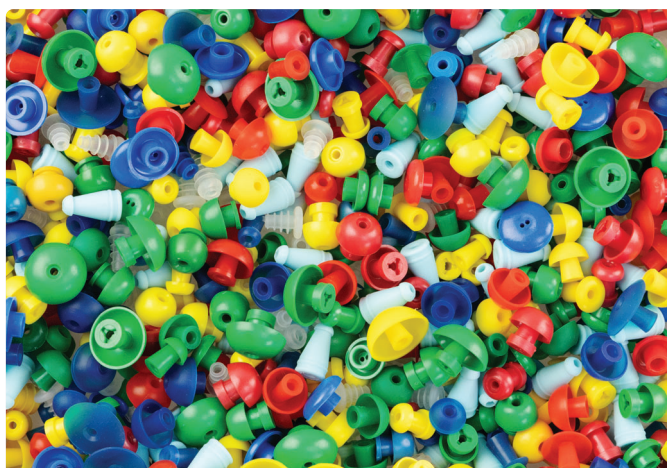

Save time pr. patient, earn more money and help more patients!

Use a fresh ear tip pr. patient and be sure to reach a seal right away, when inserting the ear tip. A fresh ear tip contains the optimal conditions for reaching a seal when inserting it into the patients' ear, whereas a multi-used one often will cause you trouble and delay you in your process. By allowing yourself to plan less time pr. patient for testing, you can see more patients pr. day. That is worth much more than cleaning an ear tip, both for your clinic and your patients.

Disposable ear tips

\$122



Otodynamics Ear Tips

Eartips for measuring DPOAE and TEOAE. For use with Otodynamics ILO-OAE instruments with DPOAE and TEOAE probes.

Part no.	Description	Unit
1160965	R7M 7mm (blue)	Bag/100 pcs
1160966	R8M 8mm (blue)	Bag/100 pcs
1160967	R9M 9mm (blue)	Bag/100 pcs
1160959	R11M 11mm (blue)	Bag/100 pcs
1160961	R13M 13mm (blue)	Bag/100 pcs

\$122



Otodynamics Ear Tips

For UGD or SND probes

Part no.	Description	Unit
1160962	T-R4.8S 4.8mm (blue)	Bag/100 pcs
1160963	T-R5.8B 5.8mm (blue)	Bag/100 pcs
1160964	T-R6.8B 6.8 mm (blue)	Bag/100 pcs

\$186



DP Probe Ear tips

Foam probe tips for paediatric and small adult ears. For use with X-series and U-series probes.

Part no.	Description	Unit
1160957	Large	Bag/100 pcs
1160958	Small	Bag/100 pcs

\$590



NEW

Clinical Audiology Kit

Suitable for students, graduates and new clinics

Part no.	Description	Unit
1161195		1 kit
Includes:		
1x Heine Mini 3000 Otoscope		
1x LED Earlight		
50x Heine Disposable Specula 4.0mm		
50x Heine Disposable Specula 2.5mm		
1x Stetoset		
1x grey syringe		
1x Combo 3 in 1 tool		
1x Vent cleaner		
5x Sanibel Universal Type A Probe tube		
20x Alcohol Pads		
20x 10mm Green impedance tips		
20x 11mm Yellow impedance tips		
20x 12mm Red impedance tips		
20x 13mm Blue impedance tips		
20x 14mm Green impedance tips		
20x 3A foam tips		
20x 3B foam tips		
20x Disposable ear cushion covers		
2x Clear plastic containers		

Probe tips

\$5



IMP

Probe tip for Sanibel™ ADI™ ear tips

Fits all Otowaves, new MT10 and GSI Allegro.

Part no.	Unit
8002592	1 pc

\$30



IMP

Probe tip for Sanibel™ ADI™ ear tips

Fits Titan™ and easyTymp™.

Part no.	Unit
8001701	1 pc



IMP OAE ABR

Probe tip kit

Fits AT235™ or newer, AA222™, Sera™, Titan vers. 2, Tympstar Pro™ Novus™, easyScreen™, and touchTymp™. Includes probe tip and gasket.

Part no.	Unit	
8107762	Set/ 2 pcs	\$33
8501570	3 sets	\$85

\$70



OAE

Internal probe tip replacements for ERO•SCAN™

Fits ERO•SCAN™.

Part no.	Unit
8102339	Bag/4 pcs

\$1414



OAE

ERI0D probe tips for ERO•SCAN™ Pro and OtoRead™

It is recommended to clean the probe tip after every patient and to change the probe tip after 20 tests (40 ears). No re-calibration necessary.

Part no.	Unit
8012836	Bag/4 pcs

\$55



OAE

Remote probe tip replacements for ERO•SCAN™

Fits ERO•Scan™ Classic.

Part no.	Unit
8102340	Bag/4 pcs

\$242



OAE

Probe tip kit for OAE Ø3.3 mm ear tips

Fits all Eclipse™ – use with ADI™ ear tips. Requires no re-calibration.

Part no.	Unit
8012960	Bag/8 pcs

\$50



IMP

Sanibel™ probe tip adaptor kit

Converts AT235™, AT235h™, AA220™, AA222™ (blue/new), MI24™, MI26™, MI34 old series, MI44™ for use with ADI™ ear tips. Includes probe tip, washer, damper. No re-calibration required.

Part no.	Unit
8500430	Set/3 pcs

Probe tips

\$60



100 pcs

OAE

Sanibel™ probe tubes for ADI™ ear tips

Fits Corti™, ERO•SCAN™ Plus and the new OtoRead™.

Part no.	Unit
8104159	Box/100 pcs

\$15



OAE

Probe Tube Removal Tool

Fits Corti™, ERO•SCAN™ Plus and the new OtoRead™.

Part no.	Unit
8513131	Bag/1 pc

\$65

NEW



1 pc +
100 pcs

OAE

Probe Tube Removal Tool kit

Fits Corti™, ERO•SCAN™ Plus and the new OtoRead™.

Part no.	Unit
8513132	Bag/100 pcs + 1 pc

Sanibel™ Probe Tube Removal Tool









A = Probe
B = Probe Tube Installer
C = Attachment to probe cable
D = Tube Remover

To Remove a Probe Tube

1. Place the front hole of the Probe Tube Tool over the end of the probe, as shown. The probe should be seated against the face of the tool and snap in place.
2. Squeeze the tool closed and hold it.
3. Twist the tool a couple of times while holding the tool closed and pulling away from the probe.
4. The probe tube will pull out from the probe. Discard the probe tube.

To Install a New Probe Tube

5. Place a new tube into the hole on the top of the tool.
6. Align the probe with the tube.
7. Push the probe onto the probe tube until the probe contacts the face of the probe tube tool and a snapping sound is heard.
8. Slide the tool off, leaving the tube in place in the probe.

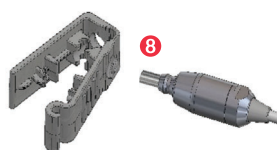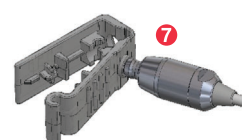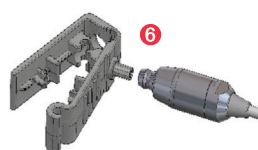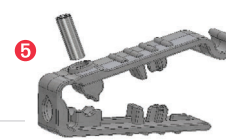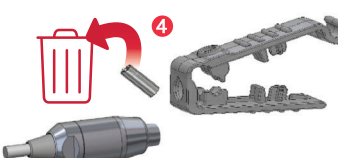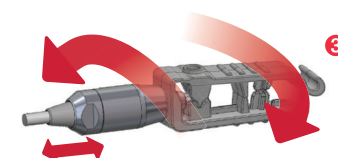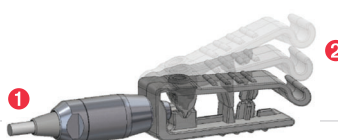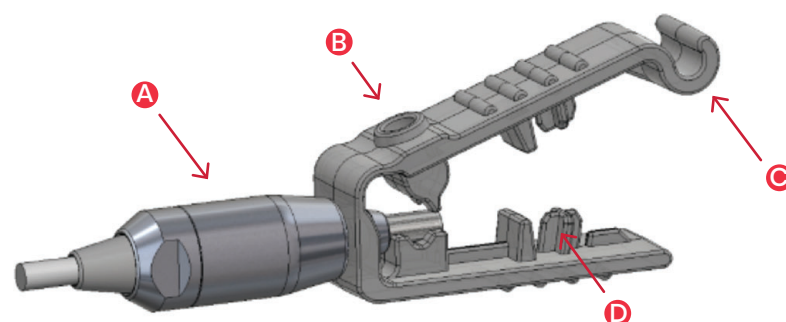

Supplies for insert earphones



ABR AUD

Sanibel™ foam tip, small

Compatible with all insert phones.
9.7mm

Part no.	Description	Unit
8527913	Pediatric pink	100 pcs
8527915	Pediatric pink	750 pcs



ABR AUD

Sanibel™ foam tip

Compatible with all insert phones.
13.7mm

Part no.	Description	Unit
8527920	Adult pink, medium	100 pcs
8527921	Adult pink, medium	750 pcs




ABR AUD

Sanibel™ foam tip

Compatible with all insert phones.

Part no.	Description	Unit
8535534	Adult pink, jumbo	24 pcs

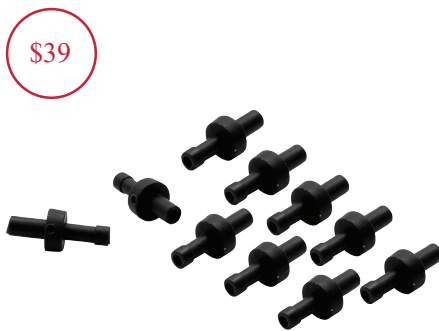

ABR

Sanibel™ IP30™ tube connectors

For use with IP30™ and E-A-Rtone® 3A.

Part no.	Unit
8107391	Bag/10 pcs




ABR AUD

Sanibel™ ADI ear tip adaptor

Makes Sanibel™ ADI ear tips compatible with 3A insert phones.

Part no.	Unit
8500580	Bag/6 pcs

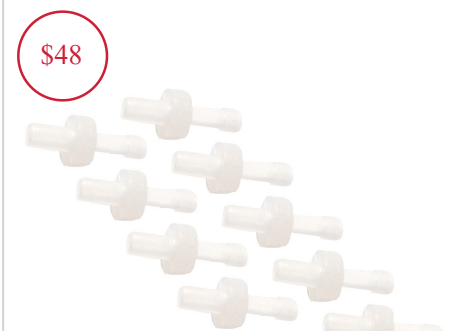

ABR AUD

Tube nipples for E-A-Rtone® 3A

For use with sound tubes 3A.

Part no.	Unit
8500140	Bag/10 pcs



ABR AUD

3A E-A-RLINK™ foam ear tips

Complies to all insert phones.
13.7mm

Part no.	Description	Unit
8106226	Adult/Medium size	Bag/100 pcs




ABR AUD

3B E-A-RLINK™ foam ear tips

Complies to all insert phones.
9.7mm

Part no.	Description	Unit
8106227	Pediatric/Small size	Bag/100 pcs

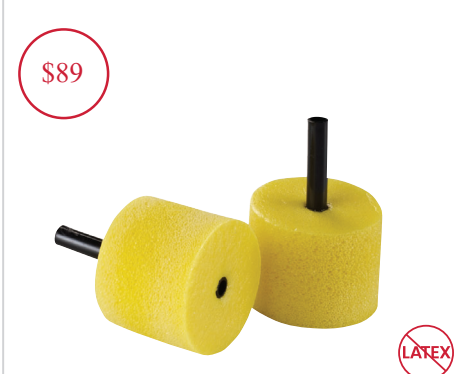

ABR AUD

3C E-A-RLINK™ foam ear tips

Complies to all insert phones.
17.8mm

Part no.	Description	Unit
8500100	Jumbo/Large size	Bag/24 pcs

Supplies for insert earphones

\$100



ABR

E-A-RLINK™ infant ear tips

For use with infants or in very small ear canals.

Part no.	Description	Unit
8500300	3.5 mm	Bag/20 pcs
8500310	4.0 mm	Bag/20 pcs

\$41




ABR

Sound tubes for E-A-Rtone® 3A 10 ohm

Includes adaptor nipples.

Part no.	Unit
8107396	Bag/1 pc

\$60

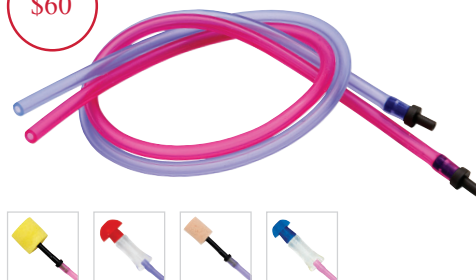

ABR

Sanibel™ IP30 and E-A-Rtone® sound tubes and nipples

Fit all IP30™ and E-A-Rtone® versions. Includes tube connectors, 1 pc blue and 1 pc red.

Part no.	Unit
8501072	Set/2 pcs

\$7



ABR OAE

Sanibel™ ear tip adaptor cleaning tool

Small brush for cleaning the Sanibel™ probe tip adaptors.

Part no.	Unit
8500830	Bag/10 pcs

\$10



HIF


Wax Brush 3 in 1

Includes:

- Full brush to remove unwanted wax from hearing aids.
- Magnet at the top of the brush to remove batteries
- Wire look at opposite end

Part no.	Unit
1160206	1 pc

\$82

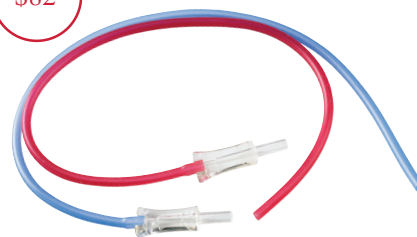

ABR AUD

Sanibel™ sound tubes and Sanibel™ ADI™ ear tip adaptor kit

Fits E-A-Rtone® 3A inserts and IP30™ insert phones. Includes ADI™ ear tip adaptors, 1 pc blue and 1 pc red.

Part no.	Unit
8100590	Set/2 pcs

\$81



ABR

Sanibel™ fastener kit for IP30™

Includes velcro neckstrap and clip.

Part no.	Unit
8507410	Set

\$13



NEW

IP30 Clip

Replacement velcro clip for insert earphones. Compatible with a variety of insert earphones, including: IP30, E-A-Rtone 3A and 5A, Etymotic ER-3A and ER-3C.

Part no.	Unit
8537295	1 pc

\$132



NEW

HIF

Foam Ear Moulds

Used as a temporary/test tip with any Post Auricular Hearing Aid (BTE). Pre-bent Clear #13 tube. Standard 13.7mm diameter.

Part no.	Unit
1161185	Bag/25 pcs



Sanibel™ Infant EarCup™ with snap or tab electrodes

Reduce test time! Test on both ears in one screening. Compatible with Titan™, AccuScreen® (1st generation), MB 11 Classic and Echo-Screen®. Use only with equipment calibrated for over-ear style disposables. Disposable over-the-ear earphones are high quality, low cost, and ideal for bilateral newborn hearing screenings. 1 set consists of 2 Infant EarCups™, 3 snap or tab electrodes and 1 alcohol pad.

For use with a range of different providers of newborn hearing screening equipment: **Grason Stadler:** Novus. **Interacoustics:** Titan, Sera. **Maico:** easyScreen, MB 11 Classic. **Madsen/Otometrics:** Accuscreen. **Natus:** Biologic Abear, Echo-Screen.

Part no.	Description	Unit	
8049600	Tab	Box/20 sets	\$396
8049700	Snap	Box/20 sets	\$396
8508848	w/o electrodes	Box/20 sets	\$302

Sanibel Infant EarCup™ w/ Hydrogel

New Hydrogel version to ensure a more comfortable solution for the baby while testing.

Choose any electrode type you wish to add to meet your preference.

Order electrodes separately.

Part no.	Unit	
8516095	Box/20 sets	\$382

NEW



ABR

Sanibel™ Infant EarCup™ adaptor set

Connects insert phones with Infant EarCups™. NB: Do not contain filters.

Part no.	Unit
8500450	Set/2 pcs

\$82

NEW



ABR

Sanibel™ sound tubes and Sanibel™ Infant EarCup™ adaptor kit

Includes adaptors, 1 blue and 1 red tube.

Part no.	Unit
8100577	Set/2 pcs



NEW

EARTurtle™ disposables ear couplers


For use with a range of different providers of newborn hearing screening equipment: **Interacoustics:** Sera, **MAICO:** easyScreen, MB 11 Classic, **GSI:** Novus, **Natus:** Echo-Screen, Echo-Screen III, Algo 3i, Algo 5, Algo 7, Biologic Abear.

No recalibration needed for any of the above-mentioned equipment.

Part no.	Description	Unit	
8529735	Tab	Box/36 sets	\$760
8530442	Snap	Box/20 sets	\$422
8522749	w/o electrodes and alcohol pads	Box/20 sets	\$353

Can be used with various equipment for ABR testing

Instruments Compatibility

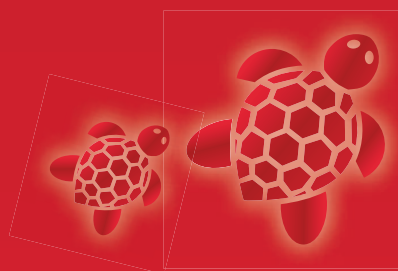

Instrument	Consumable	
	EARTurtle™	EarCup™
Interacoustics		
Titan	-	X
Sera	X	X
MAICO		
easyScreen	X	X
MB 11 Classic	X	X
GSI		
Novus	X	X
Madsen / Otometrics		
AccuScreen (1st generation)	-	X
AccuScreen (2nd generation)	-	X
Natus		
Echo-Screen	X	-
Echo-Screen III	X	-
Algo 3i	X	-
Algo 5	X	-
Algo 7	X	-
Biologic Abear	X	-
Otodynamics		
Otoport ABR	-	X
Otoport OAE+ABR	-	X
OtoNova	-	X

Sanibel Supply A/S
info@sanibelsupply.com
www.sanibelsupply.com

Sanibel
 Supply

Disposable electrodes

\$112

**Large, soft
electrode size!**



ABR

Sanibel™ snap electrodes

Pre-gelled with a very gentle hydrogel, perfect for ABR testing on infants. Size 2.2 x 3 cm. Item no. 8107137 is packed in 20 smaller pouches of 3 electrodes. Intended use areas: Newborn, diagnostic follow-up, ABR, C-Vemp.

Part no.	Unit
8107137	Box/60 pcs (20x3 pcs)



\$89

**Large, soft
electrode size!**



ABR

Sanibel™ tab electrodes

Pre-gelled cloth electrode with a very gentle hydrogel, perfect for ABR testing on infants. Size 2.2 x 3 cm. Packed in 20 smaller pouches of 3 electrodes. Intended use areas: Newborn, diagnostic follow-up, ABR, C-Vemp.

Part no.	Unit
8107410	Box/60 pcs (20x3 pcs)



\$19



Vermed Economy line 3-pack solid gel snap Electrode

Solid gel type snap electrode suitable for diagnostic follow-up. Latex free. Intended use areas: Newborn, diagnostic follow-up all ages.

Part no.	Unit
8516252	Bag/50 pcs



\$120



ABR

Gold Eartip Electrode

For ECoChG. Use with cable 8011343.

Part no.	Description	Unit
8500370	Adult 13 mm	Bag/20 pcs
8500440	Pediatric 10 mm	Bag/20 pcs



10 pcs \$499
1 pc \$71



ABR

Sanibel™ TM Electrode

For ECoChG. Use with cable 8500650.

Part no.	Unit
8011489	Box/10 pcs
8049550	Single pc




Choose the correct electrode for optimal performance

Sanibel Supply® offers various kinds of surface electrodes for ECG as well as EEG and ABR. Surface electrodes vary in size, adhesion, connectivity and frequency (test category).

The frequency plot to the above shows the electrical stimuli sent directly through the system and compared to the surface electrode (here ABR). The two lines are similar and demonstrate that the surface electrode can be used for ABR testing.

At Sanibel Supply®, we always test and validate our surface electrodes to ensure optimal performance for your testing!

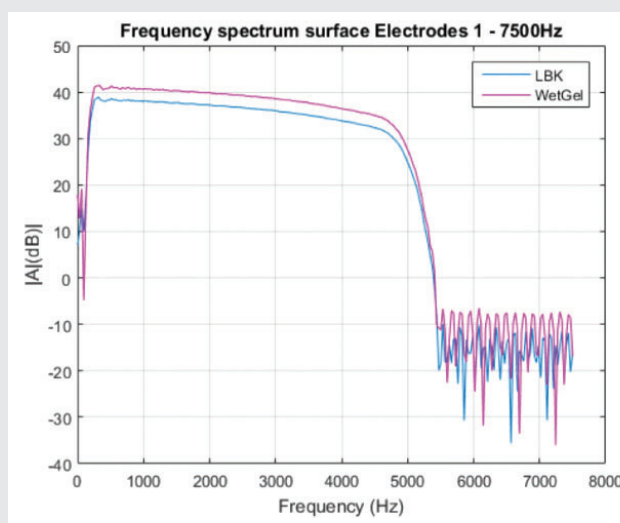

Disposable electrodes

VERMED

\$90



Vermed Wetgel Snap Electrode

Pre-gelled hydrogel electrode with Ag/AgCl sensor. Very suitable for diagnostic follow-up. Used for newborn screening, industry-wide. Packaged 3 pcs. per pouch.

Part no.	Unit
8516221	Box/75 pcs

\$45



ABR

Vermed Hydrogel Snap Electrode

Large soft electrode size making it very suitable for newborn screening and ABR measurements. Cloth material makes it convenient to wear. Size 3.4 x 2.2cm.

Part no.	Unit
8516907	Bag/60 pcs (3x20 pcs)

\$25



Vermed Tab Electrode

Rectangular tab electrode used for newborn screening, industry-wide. 12.7 x 22.9 mm. Conductive Adhesive Hydrogel. 10x10 pcs. per pouch.

Part no.	Unit
8516225	Bag/100 pcs

\$46



ABR

Vermed Hydrogel Snap Electrode

Smaller than the rectangular version, but otherwise similar. Good, soft electrode size making it very suitable for newborn screening and ABR measurements. Cloth material makes it convenient to wear.

Part no.	Unit
8102022	Bag/60 pcs (3x20 pcs)

\$22



Vermed Economy Line Snap Wet Gel Electrode

Non-woven cloth electrode with conductive wet gel. Passed ISO 10993 and is Latex free.

Part no.	Unit
8516251	Bag/50 pcs

\$17




Vermed Tab Electrode

Rectangular tab electrode used for newborn screening, industry-wide. 22.9 x 34 mm. Conductive Adhesive Hydrogel. 10x10 pcs. per pouch.

Part no.	Unit
8516227	Bag/100 pcs

AMBU

\$22

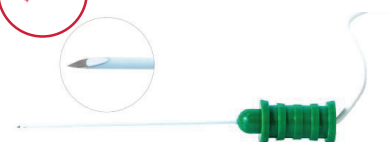

Ambu® Neurology Needle Electrode

This needle electrode features exceedingly low penetration-resistance levels. The needles have been coated with smooth coating which reduces friction when inserting and working in the target area.

Part no.	Unit
8519037	Box/40 pcs

Disposable electrodes

\$26



ABR

Ambu® 0415M tab electrode

Medium size tab electrode with solid gel for 12-lead ECG. Presented conveniently with 10 pieces to a card. Fish tail tab design reduces electrode lift-offs. Suitable for both adults and children.

Part no.	Unit
8516138	Bag/100 pcs (10x10 pcs)



\$35



ABR

Ambu® 2300 Blue Sensor Tab electrode

The medium size Ambu BlueSensor 2300 electrode features a skin-friendly solid gel that provides optimal signal quality recordings during 12-lead diagnostic testing ECG application. Very suitable for ABR diagnostics of all kinds.

Part no.	Unit
8516139	Box/1000 pcs (10x100 pcs)



\$43



ABR


Ambu® Blue Sensor™ N-00-S

Foam ECG electrode for use with button cables. Wetgel.

Part no.	Unit
8013029	Bag/25 pcs



\$8

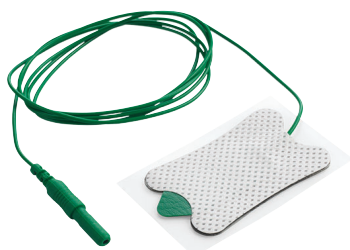

ABR

Ambu® Ground electrode 71410-M/1

Part no.	Unit
8500680	1 pc



\$31



ABR

Ambu® Neuroline 720

Ambu Neuroline 720 electrodes are pre-gelled and self adhesive. They are ideal for use in sleep registration where the strong, yet skin-friendly, adhesion ensures that the electrode stays in place for the entire examination, even if the patient sweats.

Part no.	Unit
8500060	Bag/25 pcs



\$43



ABR

Ambu® Blue Sensor™ P-00-S

High quality Ag/AgCl sensor – ensures excellent trace quality.

Part no.	Unit
8504782	Bag/50 pcs



\$50



ABR

Ambu® Neuroline 720 – with 100 cm cable and k-connector

Ambu Neuroline 720 electrodes are pre-gelled and self adhesive. They are ideal for use in sleep registration where the strong, yet skin-friendly, adhesion ensures that the electrode stays in place for the entire examination, even if the patient sweats.

Part no.	Unit
8503244	Bag/10 pcs



\$66



ABR

Ambu® Blue Sensor™ N-10-A

Foam ECG electrode with 10 cm wire and 4 mm connector, for use with universal electrode cable. Wetgel.

Part no.	Unit
8500670	Bag/25 pcs



\$10



NEW

ABR

Ambu White Sensor 40554

Specifically designed for babies and newborns.

Part no.	Unit
8520431	Bag/3 pcs





Surface electrodes tailored to your need

At Sanibel Supply® we believe that it is important to choose the right surface electrode for each specific recording and patient type.

Wetgel or hydrogel electrodes?

A wetgel electrode consists of water, thickener, ionic salts, surfactants and bactericides or fungicides. Wetgel migrates easily and rapidly onto the skin surface, reducing impedance levels, which in turn will generally result in lower noise in the recording.

A hydrogel surface electrode is less adhesive when compared to wetgel, making it particularly suitable for newborns.

Chloride salt concentration

The higher the chloride salt content, the more conductive the electrode. Higher salt content (10% and above), pre-gelled surface electrodes are useful for making fast, high quality measurements of biopotentials. High conductivity electrodes generally have reduced artifact, due to the lower impedance generated between electrode and skin surface. As the chloride salt content drops, there is a corresponding reduction in electrode conductivity.

If the chloride content is below 10%, the electrode can then be employed for longer recording times, with a reduced risk of skin irritation. As a rule of thumb, wetgel has a chloride content of 10% and hydrogel 4%.

It is important to use electrodes with similar chloride salt content when testing.

The impedance of the electrode/skin junction is highly dependent on the electrolyte type and the chloride salt concentration. For example, a hydrogel electrode with

4% chloride concentration will have approximately 10x higher impedance than a wet liquid gel electrode with 10% chloride concentration, after first application to the skin.

Adhesion

An advantage of a conductive adhesive is that the active surface can be as large as the electrode itself because no additional adhesive strip is needed. The Sanibel disposable electrodes come with a range of adhesive qualities. All electrodes are designed to adhere well to skin surfaces, but the adhesion characteristics can be grouped according to the application:

- Strong adhesive
- Moderate adhesive

Strong adhesion is not recommended for newborns.

Skin preparation

For optimum electrode to skin conductivity, the skin should be lightly brushed with a gentle abrasive gel/wipe.

An alcohol wipe is not recommended to improve conductivity as this will only dry out the skin surface.



Electrode cables

\$158



ABR

Tip trode cables

For use with TTE25 gold foil tip trodes. 1.5 mm with touch proof connectors and plated clamp clips. 50 cm long.

Part no.	Unit
8011343	Set/2 pcs

\$300



ABR

Sanibel™ ECoG starter kit

For Eclipse™ and Bio-logic® AEP system. Includes cable, gel and 2 pcs TM electrodes.

Part no.	Unit
8013209	Set/4 pcs

\$56




ABR

Sanibel™ EP electrode jumper cable

Allows two of the preamplifier inputs to be connected, e.g. for linking mastoid electrodes for cortical testing.

Part no.	Unit
8107872	1 pc

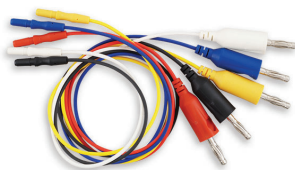

ABR

Universal electrode cable

Available in Red (R), White (W), Blue (B), Yellow (Y), Black (K). 50 cm cable. Non-tangling.

Part no.	Unit	
8011353	Set/3 pcs RWB	\$95
8011344	Set/5 pcs RWBYK	\$157

\$146



ABR

Pinch clip electrode cables

White (W), Yellow (Y), Black (K). Length 100 cm.

Part no.	Unit
8519048	Set/3 pcs WYK

\$200




ABR

Sanibel™ ECoG cable

For Eclipse™ and Bio-logic® AEP system. 1.5 mm/.060" female connector. For use with Sanibel™ TM electrodes.

Part no.	Unit
8500650	1 pc

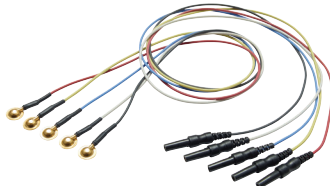


ABR

Gold-plated cup electrodes

Available in Red (R), White (W), Blue (B), Yellow (Y), Black (K). Ø 10 mm gold-plated. 50 cm cable.

Part no.	Unit	
8107418	Set/3 pcs RWB	\$212
8107145	Set/5 pcs RWBYK	\$305

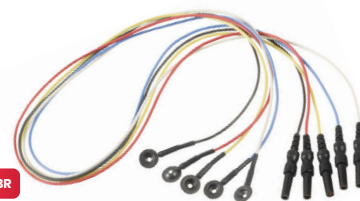

ABR

Disc electrodes

For use with Ten20™ electrode paste. Available in Red (R), White (W), Blue (B), Yellow (Y), Black (K). 50 cm cable.

Part no.	Unit	
8011498	Set/3 pcs RWB	\$186
8011345	Set/5 pcs RWBYK	\$266




ABR

Pinch style cable

For use with snap electrodes. 1.5 mm/.060" female connector. Available in Red (R), White (W), Blue (B), Yellow (Y), Black (K).

Part no.	Unit	
8500390	Set/3 pcs RWB	\$104
8500400	Set/3 pcs WYK	\$104
8500410	Set/5 pcs RWBYK	\$190

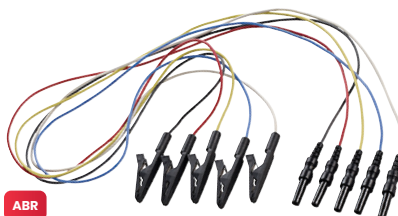

ABR

Alligator clip cable

For use with tab electrodes. 1.5 mm/.060" female connector. Available in Red (R), White (W), Blue (B), Yellow (Y), Black (K).

Part no.	Unit	
8516581	Set/3 pcs RWB	\$119
8513776	Set/3 pcs WYK	\$119
8516584	Set/5 pcs RWBYK	\$196

After use ...

After use, cables should be wiped with a soft cloth that has been dipped in a mild detergent solution or alcohol.

Contact between the cleaning solution and the connectors should be avoided.

Patient preparation

\$136



ABR

3M™ Red Dot™ trace prep, 5 m

One-step skin preparation. 18 mm x 5 m / 0.70" x 196".

Part no.	Unit
8500030	Dispenser/5 m

\$83



ABR

3M™ Transpore™ Surgical tape, 9.14 m

Transparent, latex-free, water-resistant and porous for breathability. 25 mm x 9.14 m / 1" x 10 LY.

Part no.	Unit
8500750	Box/12 rolls



\$83



ABR

3M™ Micropore™ Surgical tape, 9.14 m

Latex-free and hypoallergenic. Holds well on damp skin. 25 mm x 9.14 m / 1" x 10 LY.

Part no.	Unit
8500760	Box/12 rolls



\$20



ABR OAE

Softa swabs

70% isopropyl alcohol.

Part no.	Unit
8500050	Box/100 pcs



\$94



NEW

ABR OAE

Ambu® 2121M Skin Prep Pads

Part no.	Unit
8523414	Box/500 pcs

\$13



Gauze swabs

Non-sterile, 7.5 x 7.5 cm and 8 layered.

Part no.	Unit
8500040	Pack/100 pcs



\$13



ABR

SignaGel®, 250 g

Bacteriostatic conductive saline electrode gel.

Part no.	Unit
8102628	Tube/250 g

Patient Preparation

\$95



ABR OAE

Q-tips

Cotton tipped applicators for patient preparation. Non-sterile. Each pack contains of 100 pcs.

Part no.	Unit
8500810	Box/10 packs

2

\$17



ABR

Spectra 360® electrode gel, 250 g

Salt- & chloride-free conductive electrode gel suitable for long-term applications.

Part no.	Unit
8500740	Tube/250 g

\$73

\$131



ABR

Ten20™ conductive EEG paste, 114 g/228 g

Adhesive gel for disc electrodes. Available in pack of 3x4 oz/114 g tubes or 3x8 oz/228 g jar.

Part no.	Unit
8012924	Box/3 tubes
8500330	Box/3 jars

\$87



ABR

NuPrep™ skin prep gel, 114 g

Efficient, yet gentle skin abrasive. Chloride and acetone free. Pack of 1x4 oz/114 g tube or 3x4 oz/114 g tube.

Part no.	Unit
8500320	Box/3 tubes

\$24



ABR

Everi™ abrasive paste, 160 g

Reduces the skin impedance during EEG, EP, ECG and Holter Stress test examination.

Part no.	Unit
8012952	Bottle/160 g

\$30



ABR

Sonavelle® electrode gel, 250 ml

Formaldehyde-free and glycerol-free.

Part no.	Unit
8500380	Bottle/250 ml

Infection control

\$178

500 pcs

2

LATEX

AUD



3" Disposable ear cushion cover

Replaces disinfectants. Headphone covers for classic audiometric headphones.

Part no.	Unit
8107419	Box/500 pcs

\$60

100 pcs

2

LATEX

AUD



4.5" Disposable ear cushion cover

Replaces disinfectants. Headphone covers for noise cancelling audiometric headphones. For HDA200/DD450/DD65v2/Peltor headsets.

Part no.	Unit
8513404	Bag/100 pcs

\$20



ABR IMP OAE

Spiraflex Cleaners

Used to clean probe tips.

Part no.	Unit
1160279	Bag/30 pcs

\$7



ABR IMP OAE

Probe cleaning tool

For **Interacoustics** Titan™/AT235™/AA222™/Sera™/MT10™/**Amplivox** Otowaves/**MAICO** touchTymp MI22, MI23, MI24, MI26, MI34, MI36, MI44 + touchTymp™ and **GSI** TympStar™ Pro, Novus and Allegro.

Part no.	Unit
8011351	1 pc

\$85



Disinfectant detergent, 1000 ml

Alcohol-free, ready to use solution for a quick and extensive disinfection of non-invasive medical devices.

Part no.	Unit
8501706	Bottle, 1000 ml

\$35



Safety goggles

Part no.	Unit
1160197	1 pc

\$31



ABR IMP OAE

3-in-1 floss

Used to clean probe tips.

Part no.	Unit
8003565	Box/100 pcs

\$15



ABR IMP OAE

Bridge & Implant floss

Used to clean probe tips.

Part no.	Unit
8003566	Box/30 pcs

2

Infection control is so much more than cleaning probe tips, protecting your personal computer's keyboard, and maintaining high standards within hygiene using e.g. disinfectant detergent for a quick and extensive disinfection of non-invasive medical devices.

When it comes to infection control, preventing germs, dirt and dust from spreading via other sources, is also of critical importance – both to the safety of your patients and yourself! Gloves, shoe covers, facemasks, and aprons are just some of the "accessories", which can help you succeed in this.

Maintain high, professional standards and make your clinic an ever-safe environment for your patients, making the wearing of gloves, facemasks and other body protection wear a natural part of the "patient interaction" process.

Infection control

\$30



2

Mediwipes Antibacterial wipes canister

70% isopropyl alcohol

Part no.	Unit
1161123	160 sheets

\$35



SecurePlus Nitrile gloves

For use at medical sites and at research institutes etc. Do not interact with petroleum oils, grease, solvents. Does not contain latex. Powder free.

Part no.	Description	Unit
1161154	Extra Small	Box/200 pcs
1161177	Small	Box/200 pcs
1161167	Medium	Box/200 pcs
1161170	Large	Box/200 pcs
1161171	Extra Large	Box/200 pcs

\$14



Disposable face masks

Perform high hygiene practices in front of your patients or clients.

Part No.	Description	Unit
1160166	Safety mask with loops	Box/50pcs

\$30



PDI Detergent wipes, refill

Solo lid bag. Wipe size: 300 x 245mm.

Part no.	Unit
1160191	Bag/225 pcs

\$129



2

PDI Detergent wipes, cannister

Wipe size: 300 x 200mm

Part no.	Unit
1161196	150 pcs

\$76



KN95 Disposable facemasks

Keep a professional appearance also after servicing a hearing instrument.

Part no.	Unit
1160804	Box/80 pcs

REM probe tubes for HA fitting

\$110



DISCONTINUED
LIMITED STOCK
AVAILABLE

HIF

Sanibel™ REM probe tubes, type A
Universal fit, length 76 mm, beveled
45° angle tip and black marker ring. **LATEX**

Part no.	Unit
8104116	Tray/50 pcs 2

\$100



DISCONTINUED
LIMITED STOCK
AVAILABLE

HIF

Sanibel™ REM probe tubes, type B
Fit all brands, **Interacoustics**,
Natus Aurical, **MedRx** and **AudioScan**. **LATEX**

Part no.	Unit
8500360	Bag/40 pcs 2

\$23



50 pcs

HIF

Silicone tubes for SPL60™
For use with Affinity™ REM probe. **LATEX**

Part no.	Unit
8500700	Bag/50 pcs 2

\$84



NEW

Sanibel probe tubes, type A
Universal REM Probe Tubes

Part no.	Unit
8533301	50 pcs

\$84



NEW

Sanibel probe tubes, type C
Universal REM Probe Tubes

Part no.	Unit
8533303	50 pcs

\$84




NEW

Sanibel probe tubes, type D
Universal REM Probe Tubes

Part no.	Unit
8533304	50 pcs

\$110

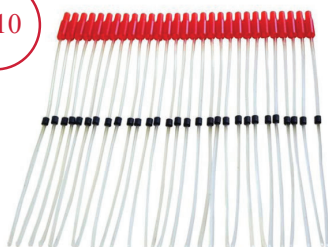

Audioscan probe tubes

Compatible with Audioscan, Axiom and
Verifit equipment

Part no.	Unit
8512115	40 pcs

Use Sanibel™
probe tubes
for any
equipment

The "REM effect"
– better fitting & fewer follow-ups!

Be sure of an
individual and
ideal fitting

Assurance of your
fitting gives you
confidence in your work

Confidence increases
your credibility and
client satisfaction

Client satisfaction means
**fewer follow-up
appointments**


We often hear from the industry that only the manufacturers' specific probe tubes should be used with their specific equipment. However, a recent comparative study by Sanibel Supply® disproves this myth.

In this study, our in-house audiologists and engineers compared the effects of different probe tubes on an Affinity²⁰ system. They compared AudioScan tubes (red), Otometrics tubes (blue), Sanibel™ tubes (orange) and Sanibel™ tubes (grey).

This involved running each tube via the 'tube calibration' process and then comparing the following REUG curve when the tube had been left in the calibration position.

The investigation showed no significant difference in the REUG curve and all probe tubes achieved the same result when the calibration process had been run correctly.

Impression Materials and supplies

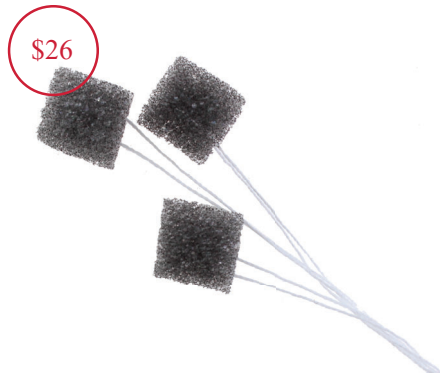

\$26

HIF

Foam Impression Pads, standard, small, grey

Diameter approx. 0.95 cm.

Part no.

1160269

Unit

Pkt/ pcs

2



\$26

HIF

Foam Impression pads, standard, medium, white

Diameter approx. 1.1 cm.

Part no.

1160862

Unit

Pkt/50 pcs

2



\$26

HIF

Foam Impression pads, standard, large, blue

Diameter approx. 1.25cm.

Part no.

1160860

Unit

Pkt/50 pcs

2



\$56

HIF

Foam Impression pads, conical, small, pink

Diameter approx. 1.1 cm, length 1.2 cm.

Part no.

1160861

Unit

Pkt/100 pcs

2



\$64

HIF

Foam Impression Pads conical standard, yellow

Cone shaped dams with thread. Comes in standard size for CIC Impressions.

Diameter approx 1.1 cm, length 0.7cm

Part no.

1160863

Unit

Pkt /100 pcs

2



\$48

HIF

Cotton Impression Pads, small, - blue

With blue safety thread.
Diameter approx, 0.65cm

Part no.

1160227

Unit

Pkt /50 pcs

2



\$48

HIF

Cotton Impression Pads, medium, white

With white safety thread.
Diameter approx, 0.75cm

Part no.

1160228

Unit

Pkt /50 pcs

2



\$48

HIF

Cotton Impression Pads, large, red

With red safety thread.
Diameter approx, 0.85cm

Part no.

1160229

Unit

Pkt /50 pcs

2

Reasons to maintain high-level skills in the art of ear impressions

The need to maintain the skills within ear impressions using quality products when working with patients is evident. Your patients will thank you for a high level of care and professionalism, long-term. Science and art come together when making a good ear impression. A deep understanding and conscious choice of the types of materials used, combined with the professional desire to be precise and demonstrate great quality and accuracy in your work, will help you make great impressions for each and every one of your patients – regardless of their different ear canal anatomy.

There is a need to be skilled at ear impressions for several other reasons;

- 1- The concept of “patient delight”: One of the most frequent components of “delight” is you – the professional. Sanibel Supply agrees with all hearing professionals viewing the skill of impression taking as a critical one that impacts the potential to delight patients.
- 2- Bad impressions result in bad earmolds or shells, resulting in bad fittings and dissatisfied patients.
- 3- The opportunity to provide custom solutions for your patient in a non-custom world is one of the key benefits of being skilled and proficient at making ear impressions.
- 4- The win-win of demonstrating your professionalism through the care and effort it takes to provide the best impression of anatomy, ear texture, and desired products.
- 5- There are a lot of other options where we might need an impression, other than when fitting non-custom technology.



Fotoplast® Gel, 20 ml

Light curing repair resin.

Part no.	Description	Unit
8530895	Clear transparent	1 pc
8530896	Opaque H (yellowish)	1 pc



Fotoplast® Lacquer 3

suitable for all hard earmolds. Curing time: 60seconds. Tissue-friendly surface sealing.

Part no.	Description	Unit
8530893	50 ml	1 pc



Detax Pink Putty Impression Kit

This impression putty has a soft creamy consistency for easy mixing and removal. The Kit includes two dosing spoons and 500ml of catalyst and 500ml of base material to produce approximately 60 impressions.

Part no.	Unit
1160519	2x500 g jars



Otoform A/K 2 part silicone

Proven kneadable impression material, non-sticky and smooth with a slightly pressure-forming consistency. An even and controlled pressure is applied onto the ear tissue. Otoform AK is advantageous for little resilient ear tissues for fitting and sealing of earmolds

Part no.	Unit
1160207	2x800 g jars
1161178	2x400g jars



HIF

Impression material A-ZOFT

Premium quality and smooth initial consistency. Extremely soft, easy to remove and allows highest impression precision. Final hardness approx. 30Shore A.

Part no.	Unit
1160282	2x800 g jars
1160283	40 foil bags (sets 2x8 g)

Impression Materials and supplies

\$10

OTO-EASE™ Earmold & Hearing Aid Lubricant

A grease-less bacteria free liquid lubricant which eases insertion of ear-molds, I.T.E and C.I.C hearing aids




Part no.	Unit
1160500	1x 15ml

\$132

2

Bite blocks

Used to stabilize the jaw for open mouth impressions.

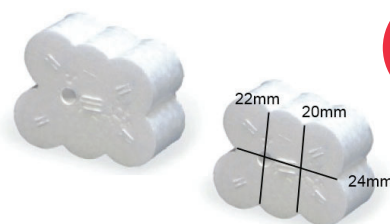

NEW

Part no.	Unit
8527849	Bag/100 pcs

\$18



Flat 4" Type Mixing Spatula

Part no.	Unit
1160192	1 pc

\$4



Mini Cleaning Brush

B & C components for ear impressions

Part no.	Unit
1160270	1 pc

\$39



Mixing tray with dosing scale

Part no.	Unit
8522105	1 pc

\$24



NEW

Mixing Pad

Disposable paper pads for mixing impression material. Foam, no slip surface on bottom of pad. Package of two pads with 35 sheets per pad.

Dimensions 6" x 6"

Part no.	Unit
8522106	Bag/ 100 sheets

\$48

\$44

\$11

Professional Impression Syringes

Easy to fill, inject and clean. Provides a stable and easy injection

Part no	Description	Unit
1160274	4mm Grey	1 pc
1160261	3mm Blue	1 pc
1160214	Economy 2lg	1 pc



Impression Materials and supplies

\$7



\$74



Earlight with Scaled Top

Part no	Description
8519357	Luxamed Earlight with acrylic scaled top, black
8519358	Luxamed acrylic top with scale

\$95

Pink Mixing Tips, Medium

Suitable with Otoform® A soft and
Otoform® A Flex
5.4mm, 50ml

Part no.	Unit
1160167	Bag/100 pcs



NEW



Surgical Blades and Handles

Suitable for cutting HA tubing.

Part no.	Description	Unit	
8515427	Surgical blade handle No. 3	1 pcs	\$29
8515428	Surgical blades No11	100 pcs	\$100

\$125



NEW

Manual Impression Gun

Injection. Used in conjunction with Otoform® AK, A flex, A soft A softX

Part no.	Description	Unit
1161029	For 50 ml double cartridges	1 pc

\$184

Blue Mixing Tips, Large

Suitable with Otoform® A soft and
Otoform® A Flex
6.5mm, 50ml

Part no.	Unit
8522183	Bag/100 pcs



NEW

\$36

Loctite 406 Instant adhesive

An instant adhesive designed especially
for the fast bonding of plastics and
elastomeric materials where very fast
fixturing is required. This product offers
exceptionally low viscosity.

Part no.	Unit
8005067	1 pc



\$17



NEW

Auragel ear mold sealant

Auragel is a safe, water-based gel suitable for use with ear moulds,
custom hearing aids and ear plugs.

Part no.	Description	Unit
1160226	20g	1 pc

Impression Materials and supplies

\$200

NEW



Otoform® A flex

Soft, flexible, easy processing. Recommended for CIC impressions, narrow auditory canals or people sensitive to pain. Final hardness approx. 30Shore A. Included 8 x S 50 double cartridges of 2 x 25 ml A+B.

Note: requires an impression gun and mixing tip.

Part no.	Description	Unit
8522433	50 ml double cartridges	1 pc

\$200

NEW



Otoform® Ak

Non-sticky, very pliable impression silicone. Pressure forming. Final hardness approx. 35 Shore A. Note: requires an impression gun and mixing tip. Included 8 x S 50 double cartridges of 2 x 25 ml A+B.

Part no.	Description	Unit
8530889	50 ml double cartridges	1 pc

\$200

NEW



Otoform® A soft

Soft, easy processing. Recommended for pressure-free impression taking and for CIC hearing systems. Final hardness approx. 40 Shore A. Note: requires an impression gun and mixing tip. Included 8 x S 50 double cartridges of 2 x 25 ml A+B.

Part no.	Description	Unit
8536593	50 ml double cartridges	1 pc

\$200

NEW




Otoform® A softX

Soft, easy processing with optimal fluidity. Recommended for deep ear impressions and sensitive ears. Final hardness approx. 25 Shore A. Note: requires an impression gun and mixing tip. Included 8 x S 50 double cartridges of 2 x 25 ml A+B.

Part no.	Description	Unit
8536594	50 ml double cartridges	1 pc

Hearing aid fitting and supplies

\$25

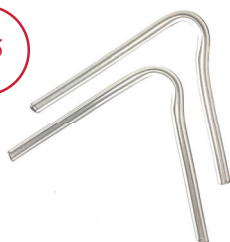

Pre-bent Tubing

Quilled on one end for easy insertion into the earpiece and double bent to lie close to the head.

Part no.	Description	Unit
1160109	Standard #13	Bag/10 pcs
1160110	Medium #13	Bag/10 pcs
1160111	Thick #13	Bag/10 pcs
1160503	Extra Thick #13	Bag/25 pcs

\$24



Moisture free Pre-bent tubing

Part no.	Description	Unit
1160122	Standard #13	Bag/10 pcs
1160123	Double wall #13	Bag/10 pcs

\$83



NEW

Tubing Retention System (TRS)

Designed for all soft earpiece materials. Each TRS tube uses an acrylic ring to hold it in position in the earpiece. TRS tubes are easily replaced and will not tear silicone materials.

Part no.	Description	Unit
8527927	Thick #13	Bag/25 pcs
8527928	Super Thick #13	Bag/25 pcs

\$25



Pre-bent Tubing with Tube Lock

Tubing system that uses a brass retention ring for soft earpieces.

Part no.	Description	Unit
1160112	Thick #13	Bag/10 pcs

\$42



Moisture Free Libby Horn

Special tube designed to increase gain over conventional tubing. Ideal when using an analog hearing instrument and fitting a high frequency hearing loss. Also designed to allow for a friction fit into earmold without the need to damage soundbore when removing it.

Part no.	Description	Unit
1160120	3 mm	Bag/10 pcs
1160121	4 mm	Bag/10 pcs

\$37



Libby Horn

Special tube designed to increase gain over conventional tubing. Ideal when using an analog hearing instrument and fitting a high frequency hearing loss. Also designed to allow for a friction fit into earmold without the need to damage soundbore when removing it.

Part no.	Description	Unit
1160118	3 mm	Bag/10 pcs
1160119	4 mm	Bag/10 pcs

\$44



NEW

Pre-bent tubing

Part no.	Description	Unit
8530962	Ø 2.0 x 3.0mm	Bag 10 pc

\$3



Tube Lock without tube

Brass tube-lock for use with 13 thick and 13 super thick tubing.

Part no.	Description	Unit
8527858	Thick #13	1 pc
8527859	Super Thick #13	1 pc

\$37



Pre-bent Tubing with Tube Lock

Tubing system that uses a plastic retention ring for soft earpieces.

Part no.	Description	Unit
1160124	Thick #13	Bag/10 pcs

Hearing aid fitting and supplies

\$37



Pre-bent Tubing with Tube Lock

Tubing system that uses a **small** plastic retention ring for soft earpieces.

Part no.	Description	Unit
1160125	Thick #13	Bag/10 pcs

\$49



Pre-bent Tubing with Tube Lock

Moisture free tubing system that uses a plastic retention ring for soft earpieces.

Part no.	Description	Unit
1160126	Thick #13	Bag/10 pcs
1160127	Standard #13	Bag/10 pcs

\$49



Pre-bent Tubing with Tube Lock

Moisture free Tubing system that uses a **small** plastic retention ring for soft earpieces.

Part no.	Description	Unit
1160128	Thick #13	Bag/10 pcs

\$63

NEW



Arc Tubing connection with retentions

Part no.	Unit
8530938	Bag 10 pcs

\$12



Tubing Puller

For replacing worn out tubing. Pulls tubes easily through any type of mold or canal tip

Part no.	Unit
1160208	1 pc

\$32

NEW



OtoVita Tubing

Designed for all soft earpiece materials. Each TRS tube uses an acrylic ring to hold it in position in the earpiece. TRS tubes are easily replaced and will not tear silicone materials.

Part no.	Description	Unit
8530957	Floss micro	30 pcs
8530958	Floss	30 pcs

\$44

NEW



Tubing

Part no.	Description	Unit
8530939	Ø 2.0 x 3.0mm	Roll 10 m

\$58

\$9

\$15



Stethoclip

Hearing Instrument listening stethoscope for reassuring the fitting process

Part no.	Description	Unit
8515600	Stethoclip, black	1 unit
8520346	Ear funnels	1 unit
8520347	ITE Adaptor	1 unit

\$46

NEW



Lightweight plastic stethoscope


For use in the testing of hearing aids. The test indicates whether the aid is weak and/or distorted. Body aids may also be tested by detaching the tubing and snapping the receiver directly into the socket. Comes with tubing and snap-in end. Ear tips are removable for cleaning or replacement.

Part no.	Description	Unit
8515600	clear, 26.2g	1 unit

Clinical Tools

\$37

NEW

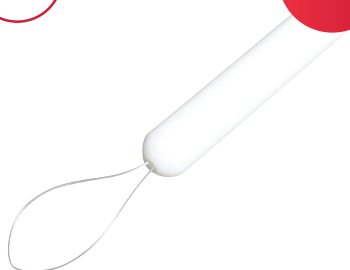

Tube Puller, Monofilament

Soft and flexible enough to easily slip through the sound bore of an acrylic earpiece, allowing quick and simple tube replacement.

Part no.	Unit
8527877	1 pc

\$37

NEW



Tube Puller, Wire Loop

Rigid enough to slip easily through the sound bore of a vinyl or silicone earpiece, allowing quick and simple tube replacement.

Part no.	Unit
8527878	1 pc

\$37

NEW



Reamer for 13 Tube

Part no.	Unit
8536590	1 pc

\$88



Probe Set, Stainless steel

For probing, material removal.

Part no.	Unit
8527929	Set of 4 pcs

\$18



Scissors, 3.5"

Sharp scissors suitable for general use in the clinic.

Part no.	Unit
8527905	1 pc

\$7



Mini Pliers

Part no.	Unit
1160222	1 pc

\$7



Mini wire cutters

Part no.	Unit
1160223	1 pc

\$10



Cleaning tool

Part no.	Unit
1160280	Bag/20 pcs

\$47

NEW



Mirror Go Round

Part no.	Unit
8536591	1 pc

Clinical Tools

\$106



Tubing Expander Plier

Stainless steel pliers with two-way lever spring allow you to quickly and easily stretch acoustic tubing before fitting it over earhooks, elbows, etc.

Part no.	Description	Unit
8527876	Short Straight	1 pc

\$135



Alligator Foreceps

Standard 5.5" stainless steel alligator forceps with serrated jaw. Not intended for invasive medical procedures.

Part no.	Unit
8527901	1 pc

\$153

Price drop



Nelson tool

Part no.	Unit
1161172	Set



NEW

Tweezers

For a variety of clinical uses. Good for working on hearing instruments and other small parts.

Part no.	Description	Unit	
8527902	Curved Point Precision Tweezers 6"	1 unit	\$27
8527903	Fine Tipped 4.75" Tweezers	1 unit	\$34
8527904	Curved Point Precision Tweezers 6"	1 unit	\$35
8536603	Bayonet Tweezers	1 unit	\$56

\$19



\$27



NEW

Precision Screwdriver set

Screwdriver set for making adjustments to hearing instruments or other precision work.

Part no.	Description	Unit
8527896	Four flat-tip screwdrivers and two Phillips tips.	Set/6 pcs
8527897	Six flat-tip screwdrivers, three Phillips tips, one needle-point tip, and one magnet tip.	Set/11 pcs

\$31



Sprue Cutter, 4.5"

A unique design easily trims tubing flush to the earpiece.

Part no.	Unit
8527935	1 pc

\$31



Needle Nose Pliers 5"

Needle nose pliers for a multitude of uses in the hearing healthcare clinic.

Part no.	Unit
8527898	1 pc

\$31



End Nipper Pliers

A perfect tool for trimming tubing from the end of a canal after the earmold has been retubed.

Part no.	Unit
8527900	1 pc

Hearing aid maintenance and hearing care

Take care of your hearing aids

Take good care of your hearing aid. Your Hearing instrument is a high tech device needing regular care in order to provide you the best performance every day.

\$7

NEW



Cleansing tablets

Part no.	Unit
8520148	Pk/4 pcs

\$9

NEW



Dehydrating capsules

Part no.	Unit
8520149	4 pcs

\$4

NEW



Drying cup

Part no.	Unit
8520150	1 pc

\$4

NEW



Cleaning cup

Part no.	Unit
8520151	1 pc

\$29

NEW



Otofloss

Part no.	Unit
8520152	1 pc

\$131

NEW



Perfect Dry LUX box for all types of hearing aids

Part no.	Unit
8520154	1 pc

\$28

NEW



Starter kit

Part no.	Unit
8520155	1 kit

\$15

NEW



Drying spray

Part no.	Unit
8520153	1 pc

\$7

NEW



Earmold Air Blowers

Part no.	Unit
1160501	1 pc

Hearing aid maintenance and hearing care

\$8

NEW



Cleaning spray

Part no.	Unit
8520147	1 pc

\$10

NEW



Cleaning wipes

Part no.	Unit
8520146	Pk/30 pcs

\$24

NEW



FlightSafe Earplugs Fly

Sanohra Fly ear plugs. Reduces ear pressure experiences when flying.

Part no.	Description	Unit
8519351	EarPlugs Fly Kids, EN	1 set
8519352	EarPlugs Fly Adult, EN	1 set

\$7

NEW



Protect your ears

Part no.	Description	Unit
8520053	Dermoprotective Gel	1 pc

\$19

NEW



Protect your ears

Part no.	Description	Unit
8520063	Earspray	1 pc

\$19

NEW



SwimGuard Earplugs Swim

Sanohra Swim earplugs. Enables people to swim without getting water in the ear, reducing ear infections.

Part no.	Description	Unit
8519353	EarPlugs Swim Kids, EN	1 set
8520143	EarPlugs Swim Adult, EN	1 set

\$8

NEW



Protect your ears

Part no.	Description	Unit
8520052	Natural lubricant	1 pc

Cerumen Management accessories

\$150



Bionix® Safe Ear Currettes

Bionix® Safe Ear Currettes allow healthcare providers to safely and gently remove cerumen with styles to match different types of wax and personal preference. Bionix® Safe Ear Currettes are flexible, to reduce injury to the ear canal, making them safer than stainless steel currettes. Seven color coded tip styles to meet specific needs, giving you “the right tool for any job” for every cerumen situation. Also available as a variety pack.

Part no.	Description	Unit
8515585	Green MicroLoop®	50 pcs/box
8515586	Orange ControLoop®	50 pcs/box
8515587	Red AngleLoop®	50 pcs/box
8519359	Purple VersaLoop®	50 pcs/box
8519360	Blue InfantScoop®	50 pcs/box
8519372	Yellow CeraSpoon®	50 pcs/box
8519361	White FlexLoop®	50 pcs/box
8519362	Variety Kit, (10 each style & 15 white)	

2

\$210



Bionix® Lighted Ear Currettes

Bionix® Lighted Ear Currettes combine illumination and magnification to help healthcare providers take less time to see and safely remove cerumen from the ear canal. The Lighted Ear Curette with magnification brilliantly illuminates the ear canal, eliminating the practice of blindly curettaging the sensitive area. Complete visualization shortens procedure time and increases efficacy. Box of 50 pcs.

All versions include 1 Light Source and 1 magnifying lens.

Part no.	Description
8519364	Lighted AngleLoop®, 4mm diameter
8519365	Lighted MicroLoop®, 3mm diameter
8519366	Lighted WaveCurette™, 2.5mm wide tip
8519367	Lighted FlexLoop®, 4mm wide tip
8519368	Lighted CeraSpoon®, 4mm wide tip
8519369	Lighted InfantScoop®, 2mm wide tip
8519373	Lighted VersaLoop®, 3mm wide tip
8519370	Lighted CeraPik™, 2.5mm concave tip. Great for hard, impacted cerumen and foreign bodies.

2

Why does ear wax accumulate?

Disruption of the self-cleaning mechanism of the ear causes cerumen to accumulate, causing different problems within our ears. These can be related to e.g. general health of the ear as well as hearing.

Anatomical variations, infectious or dermatological diseases can cause narrowing or obstruction of the ear canal. This can then interfere with the normal migratory process of foreign agents from the ear.

Indications for treatment and treatment options

When symptomatic, cerumen accumulation requires treatment. Usually, ear wax is defined as impacted, if it disrupts visualization of the tympanic membrane, as well as being connected with symptoms.

Often, removal of the cerumen is also done to allow a needed view of the tympanic membrane, e.g. for diagnostic purposes and audiometry.

The first-line option for treatment of

symptomatic cerumen in patients is cerumenolytics (softening agents, e.g. eardrops). The use of cerumenolytic agents, such as eardrops, increase the probability of cerumen clearance.

Patients with heavier impaction or ear canal disease may require irrigation or manual removal under microscope by a trained doctor or an ENT specialist.

Cerumen Management accessories

\$210



Lighted Ear Curette Variety Pack

Part no.	Description
8519363	With 1 light source and magnifying lens. Includes 8 pcs. of the Lighted Flexloop and 7 pcs. of each of the other styles (not including Cerapik).

\$230



Lighted Forceps

Bionix Lighted Forceps combine illumination and magnification to help healthcare providers take less time to see and safely remove ear canal and nasal passage obstructions.

Part no.	Description
8519392	10 forceps, 1 light and 1 magnifying lens

\$230



Lighted articulating curettes

The lighted articulated curette combines illumination, magnification and a 90° articulating tip.

Part no.	Description
1160091	Box/50 pcs

\$420



Lighted Curette Handle

Provides the ability to quickly remove cerumen or foreign bodies. Thumb controlled suction release on handle provides procedural safety.

Box of 25 suction handles, 1 light source and 1 magnification lens.


Part no.	Description
1160092	1 light source



Cerumen Removal Tool

Part no.	Description	Unit	
8527907	w/Brush and magnet	1 pc	\$7
8527908	w/Brush	1 pc	\$5
8527909	tool	1 pc	\$4

\$5

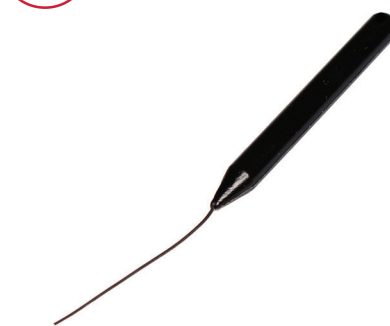

Vent Cleaning Tool

Part no.	Unit
8527910	1 pc

\$30



Tech-Care LED Ear light

The 18 Lumens light is concentrated in a small tip that fits most ears. Includes 2 AAA batteries.

Part no.	Description	Unit
1160225	Ear light with batteries and tip	1 pc
1161183	Tech-Care ear light tip with slotted base	1 pc

\$5



LED Earlite tips

To be used with LED Earlite (8527868)

Part no.	Description	Unit	
8527862	Straight tip	Bag/4 pcs	\$15
8527863	Pros tip	Bag/4 pcs	\$30
8527864	Scoop tip	Bag/4 pcs	\$22
8527865	Steel loop tip (S/L)	Bag/2 pcs	\$54

Cerumen management accessories

\$212




2


EarWay Pro – Wax Removal Helix

Very easy and safe way to remove ear wax. The helix grabs around the ear wax lump and pulls it out in a comfortable way. Available in different languages, please contact us for details.


Part no.	Description	Unit
8518025	Pro M, 5,4 mm, purple	25 units
8518024	Pro S, 4,2 mm, blue	25 units
8521854	Pro L, 6,5 mm, green	25 units

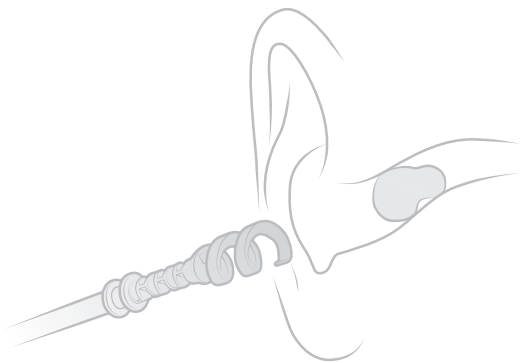


Before insertion

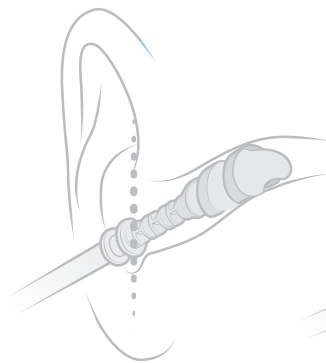

Engulfing of Cerumen

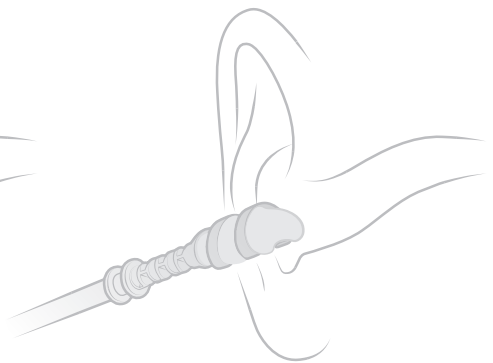

Extraction

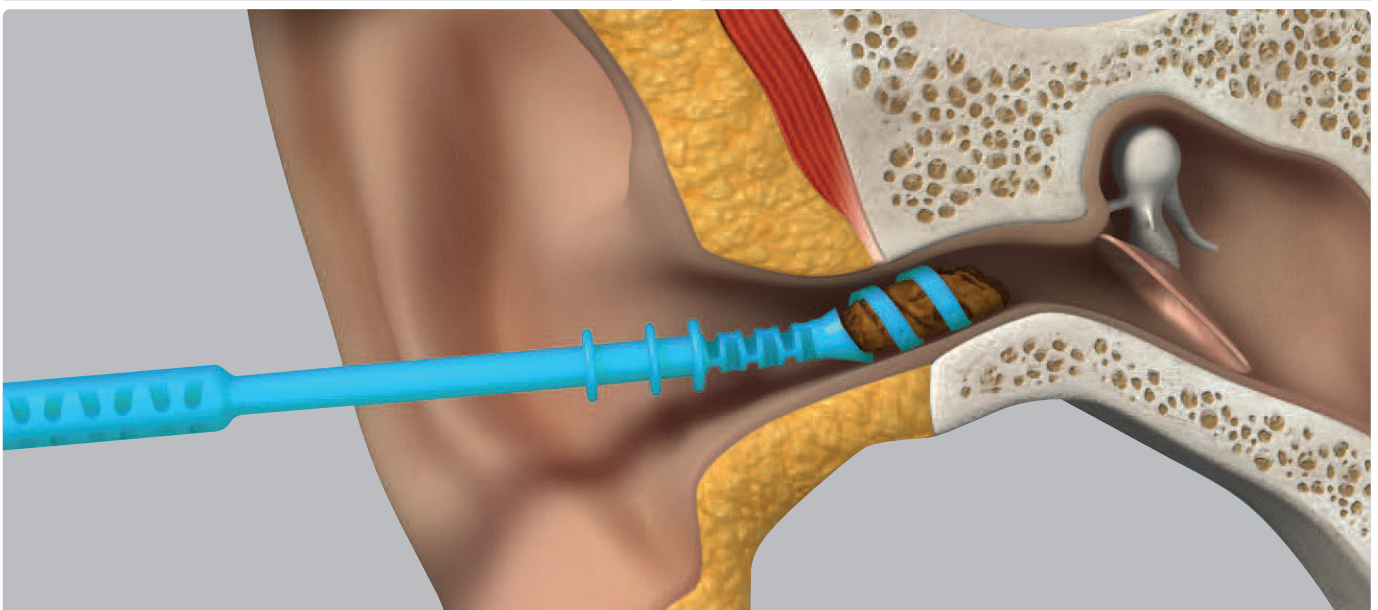

Cerumen management accessories

\$1,500



FADING PRODUCT

Minston Headlight and Loupe

Minstons Headlight and Loupe combine to give an ultra lightweight and comfortable diagnostic experience. The adjustable magnification of 2.5 - 3.5 x makes it suitable for a wide range of applications including ENT, general examination and minor procedures

Part no.	Unit
1160213	Set

Features

- * Ergonomic Design
- * Variable magnification 2.5-3.5x
- * 10,000- 40,000 Lux
- * Patented Heat Dispensing Technology
- * L Shaped Lamphead design
- * Ultra Light & Free
- * Even-light
- * Brightness control
- * Dual Battery
- * Easily adjustable

\$10



NEW

Tech care Ear wax removal kit

Complete system used for the softening, loosening, and removal of excess and unwanted earwax. The kit includes earwax removal drops and a soft rubber bulb syringe

Part no.	Unit
1160161	1 pc

\$16



NEW

Earclear ear wax remover, 12ml

Ear Clear Ear Wax Remover aids the easy removal of ear wax, helping to relieve ear wax discomfort.

Active Ingredient: Carbamide peroxide 65mg/mL

Non active ingredients (excipients): Glycerol, Succinic acid

Part no.	Unit
1160263	1 pc

Cerumen management accessories



\$1068

Tongye MC600A Suction Unit 240V

The MC600A Suction Unit is a high quality AC powered suction unit featuring step less suction adjustment, overflow protection and anti-bacterial filter. This unit is suitable for many applications such as aural micro suction, resuscitation and home care.

Part no.	Description
1160133	MC600A Suction Unit

\$16



Tongye Disposable Suction Unit Liners

To suit MC600 units. Provides a disposable and cleaning free solution for the collection of many different matter and is situated inside the hard reusable cannister included with the Tongye MC600 series suction units.

Part no.	Description
1160248	Disposable suction unit liners

\$8



Suction Tube 3M

A 3 meter silicone tube to suit the MC600 series suction units.

Part no.	Description
1160249	5mm x 3m

\$124



Tongye Overflow Stopper

Overflow stopper for Tongye Suction Units to fit all MC600 Tongye suction units

Part no.	Description
1160250	Overflow stopper

\$186



Tongye Reusable Suction Cannister

Provides collection of many different matter and is suitable for use with the MC600 series suction units

Part no.	Description
1160251	Suction cannister

Cerumen management accessories

\$7



Tongye Hose Barb

Part no.	Description
1160252	Hose Barb for MC600

\$124



Suction Controller Sterile

A small piece that goes in line between a suction tube and silicone tube that allows control of suction pressure at the fingertip

Part no.	Unit
1160253	Bag/10 pcs

\$76



Tongye Filter

Part no.	Description
1160254	Filter for MC600

\$207



Ear Speculum

Speculums have an Oval shaped tip to match the shape of the ear canal, with a thin 0.5mm wall thickness. Light grey in colour which matches the inside stainless steel speculum and is the best balance between light reflection and contrast.

Part no.	Description	Unit
1161073	Extra Small, diameter 3 x 3.5mm	Bag/100 pcs
1161173	Small, diameter 4 x 4.5mm	Bag/100 pcs
1161074	Medium, diameter 5 x 5.5mm	Bag/100 pcs
1160140	Large, diameter 6 x 6.5mm	Bag/100 pcs
1161075	Extra Large, diameter 7.5 x 8mm	Bag/100 pcs

\$124




Ear Speculum Gruber

Speculums have an Oval shaped tip to match the shape of the ear canal. Chrome plated

Part no.	Diameter - Height	Unit
1160141	25mm - 36mm	Bag/4 pcs
	26mm - 39mm	
	25mm - 39mm	
	26mm - 38mm	

Disposable Suction Tube

A single use disposable tube for the suction of fluid and other media from the ear

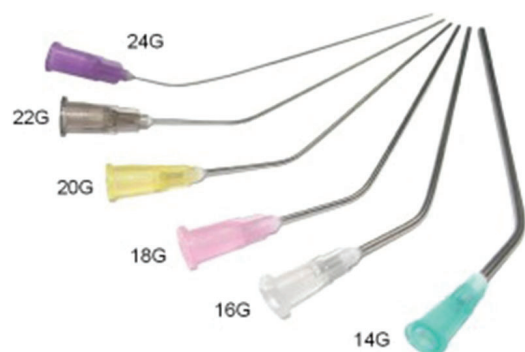

Part no.	Description	Unit
1160255	14 Gauge Bent, non-sterile Outside diameter 2.1mm	Bag/100 pcs \$225
1160256	16 Gauge Bent, non-sterile Outside diameter 1.6mm	Bag/100 pcs \$225
1160257	18 Gauge Bent, non-sterile Outside diameter 1.2mm	Bag/100 pcs \$225
1160212	20 Gauge Bent, non-sterile Outside diameter 0.9mm	Bag/100 pcs \$225
1161071	22 Gauge Bent, non-sterile Outside diameter 0.7mm	Bag/100 pcs \$407
1161072	24 Gauge Bent, non-sterile Outside diameter 0.5mm	Bag/100 pcs \$407

Balance supplies

\$108

NEW



BAL

Disposable VNG goggle foam pads

For **Interacoustics**® VNG mask.

Part no.	Unit
8011527	Box/24 pcs

2



\$17




Disposable ear specula

For AirFx™.

Part no.	Description	Unit
8013013	2.75 mm	Bag/34 pcs
8013014	4.25 mm	Bag/34 pcs

2

\$100

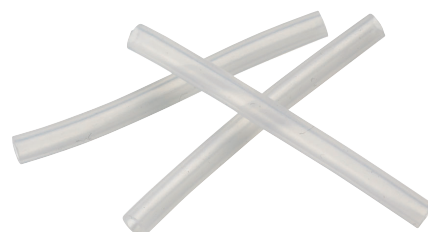

Water irrigator tips

For use with **Interacoustics**® Aqua Stim™.

Part no.	Unit
8013018	Bag/50 pcs

2

Paper

\$10



Audioscan printer paper

For VFI, Axiom and RM500SL
Width: 80 mm

Part no.	Unit
8521224	1 roll

\$10



IMP OAE

Thermal paper for Sanibel™ MPT-II printer, 20 m

For the Sanibel™ MPT-II printer which is supplied with **MAICO®** ERO-SCAN™, **Grason-Stadler®** Corti™ and **Interacoustics®** OtoRead™.

Width: 57 mm
Inner diameter: 9 mm
Outer diameter: 40 mm

Part no.	Unit
8029305	1 roll/20 m

\$20



GSI Thermal paper

For GSI39
Width: 111 x 51 x 12 mm

Part no.	Unit
8012939	1 roll

\$10



IMP

TPR10™ thermal paper for MTP10™, 25 m

Width: 79 mm
Inner diameter: 12 mm
Outer diameter: 45 mm

Part no.	Unit
8500220	1 roll/25 m

\$10



AUD IMP

TPR26™ thermal paper for Interacoustics®, 25 m

Fits most popular audiometers and tympanometers.

Width: 112 mm
Inner diameter: 12 mm
Outer diameter: 45 mm

Part no.	Unit
8500210	1 roll/25 m

\$10



IMP

Thermal paper LPQ80, 80 m

For **GSI** TympStar™ Pro.
Width: 80 mm
Inner diameter: 13 mm
Outer diameter: 80 mm

Part no.	Unit
8509893	1 roll/80 m

\$10



ABR IMP OAE

Sanibel™ MPT-III thermal paper, 20 m

For Sanibel™ MPT-III printer which is supplied with **Interacoustics®** AT235™.

Width: 79 mm
Inner diameter: 12 mm
Outer diameter: 45 mm

Part no.	Unit
8500610	1 roll/20 m

\$10



OAE

Thermal paper for ERO-SCAN™/OtoRead™, 25 m

Width: 57 mm
Inner diameter: 12 mm
Outer diameter: 45 mm

Part no.	Unit
8500200	1 roll/25 m

\$10



For OtoRead™ PRT MPT-II printer, 10 m

Width: 57 mm
Inner diameter: 8 mm
Outer diameter: 30 mm

Part no.	Unit
8013009	1 roll/10 m

Paper

\$10



ABR IMP OAE

API300 thermal paper for Titan™, 10 m

Width: 57 mm
Inner diameter: 12 mm
Outer diameter: 35 mm

Part no.	Unit
8500230	1 roll/10 m

\$10

Label paper roll HM-E200

Label size:
Width: 56 mm,
Length: 60 mm
R angle: 2



Roller size:
Paper width: 57 mm
Label spacing: 6 mm
Core diameter: 13 mm
Roller diameter: less or equal to 40 mm

Part no.	Unit
8519836	1 roll/120 labels

\$12




Audiogram chart AF12

100 pcs

Part no.	Unit
8012858	1 pad

\$13

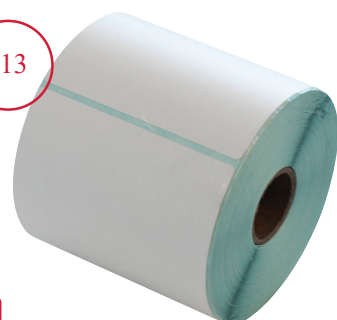

IMP

Adhesive thermal paper LPQ80, 195 labels

For **GSI** TympStar™ Pro.
Width: 80 mm
Inner diameter: 12 mm
Outer diameter: 80 mm

Part no.	Unit
8509893	1 roll/195 labels

\$25




IMP

Adhesive thermal paper LPQ80, 195 labels

For **GSI** TympStar™ Pro.
Width: 82 mm
Inner diameter: 25 mm
Outer diameter: 88 mm

Part no.	Unit
8509894	1 roll/195 labels

\$18

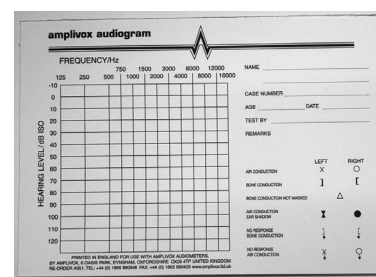

Amplivox audiogram chart

50 pcs

Part no.	Unit
8013007	1 pad

\$16




GSI Thermal paper

For GSI70
Width: 109.5mm

Part no.	Unit
8500240	1 roll

\$66

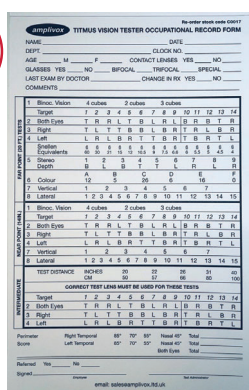


Titmus Vision Tester forms

50 pcs

Part no.	Unit
8507192	1 pad

\$66

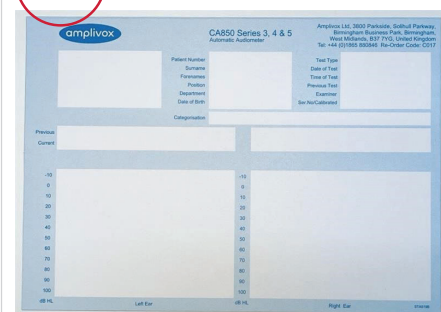

CA850/3/4/5 paper Amplivox

250 pcs

Part no.	Unit
8507761	1 pad

Otoscopy supplies

\$105



AUD IMP

Specula wall dispenser

Stylish wall dispenser for ear funnels. Can contain 6 different sizes.

Part no.	Description	Unit
8515678	Specula wall dispenser	1 pc

\$25



AUD BAL IMP

Video Otoscope specula

KaWe, fits the new **Interacoustics**® Viot and OriVision otoscope.

Part no.	Description	Unit
8515083	2.5 mm Speculum	Bag/100 pcs
8511728	4.0 mm Speculum	Bag/100 pcs

2

34 pcs \$6
850 pcs \$124

No sharp edges



AUD IMP

Welch Allyn otoscope specula

Disposable ear speculums for **Welch Allyn** otoscopes. Suitable for MacroView and Diagnostic

Part no.	Description	Unit
1160220	2.75 mm Speculum	Bag/34 pcs
1160219	4.25 mm Speculum	Bag/34 pcs
1161138	2.75 mm Speculum	Bag/850 pcs
1161122	4.25 mm Speculum	Bag/850 pcs

2

\$40



Otoscope bulb

Replacement bulb for Heine mini3000 Otoscope. 2.5v

Part no.	Description	Unit
1160276	Halogen bulb	1 pc

\$32



AUD BAL IMP

Disposable speculas for video otoscope

Fits the AURICAL™ OTOcam 300, **Interacoustics**® VIOT and XION™ video otoscopes.

Part no.	Unit
8500470	Bag/50 pcs

2

\$25

No sharp edges



AUD IMP

Welch Allyn otoscope specula

Disposable ear speculums for **Heine** otoscopes.

Part no.	Description	Unit
1160273	2.5 mm Speculum	Bag/50 pcs
1160272	4 mm Speculum	Bag/50 pcs

2

PRICE DROP

\$13



AUD IMP


Otoscope specula

Disposable ear speculums for **Heine** otoscopes

Part no.	Description	Unit
8514368	2.5 mm Speculum	Bag/100 pcs
8514367	4.0 mm Speculum	Bag/100 pcs

2

\$23

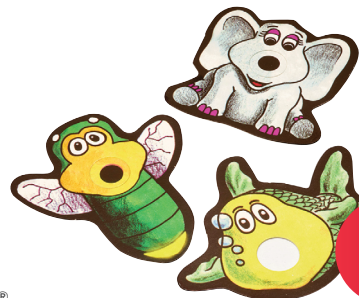

5 pcs

EARnimals®

Cartoon faces that are placed on the otoscope. Make ear examinations more fun for kids.

Part no.	Unit
8011501	Bag/5 pcs (assorted)

Otoscopy supplies

\$220



AUD IMP

HEINE® Mini3000® Otoscope

Classic pocket otoscope from Heine in Germany.

Part no.	Unit
1100218	1 pc

AUD IMP

Luxamed Otoscope Auris LED 2.5V

Otoscope from German Luxamed. LED light to ensure longevity and great light source, dimable. Fits Heine/KaWe compatible speculums.

Part no.	Unit
8519433	1 pc

\$230



Luxamed Otoscope Auris LED 3.7V USB/PC Charged

Otoscope from German Luxamed. LED light to ensure longevity and great light source. USB charging. Fits Heine/KaWe compatible speculums.

Part no.	Unit
8519434	1 pc

\$340

Accessories

\$65

NEW



Demo Ear skin coloured

Part no.	Description	Unit
8536602	Right	1 pc
8536601	Left	1 pc

\$75



Welch Allyn 3.5V Lamp

For Welch Allyn Macro-View otoscope.

Part no.	Unit
1160258	1 pc

\$114



Aidapter for hearing aid fitting

Clear Silicone Acoustics couplers. A unique silicone piece that clicks onto the HA-1 adapter.

You can use the Aidapter with receiver-in-canal (RIC) and thin tube hearing aids.

Compared to acoustic putty, the Aidapter™ provides a good acoustic seal and standardizes the coupling, ensuring repeatable coupler measurements.

The Aidapter™ also eliminates the hassle and hygienic issues of using acoustic putty.

Part no.	Description	Unit
8525568	Interacoustics	Bag/3pcs
8525578	MedRx	Bag/3pcs



We know that, first and foremost, great quality
is at the top of your list, and you'll find it in each
and every purchase
you make from Sanibel Supply.

Bringing you the best is truly what drives us day
after day.

If we can ever be of assistance, don't hesitate to
let us know, info@diatec-diagnostics.co.nz

To order, please send us an e-mail
orders@diatec-diagnostics.co.nz
or contact us on 02 8899 1277.

Sanibel Supply is a company within the William Demant Group

Sanibel
Supply

Disclaimer

This catalogue has been published for Australia wide distribution. Some parts described may NOT BE AVAILABLE in your area.
For advice on compatibility and instructions of use, please consult your operation manual or local distributor. Sanibel Supply AUS assumes no legal liability or responsibility for the accuracy, typing errors or mistakes, correctness or completeness of any information in this flyer. The information in this flyer was correct to the best of our knowledge at the time of printing. Sanibel Supply AUS reserves the right to make changes at any time in specification of products described in this flyer without notice or obligation.
MADSEN Acouscreen® and MADSEN OTOflex® are registered trademarks of Otometrics A/S. Echo-screen®, Audix®, Navigator® Pro, and ABAse® are registered trademarks of Natus Medical Incorporated. Sanibel® is a registered trademark of Sanibel Supply®. All other instrument brands or model names are copyrighted or trademarked by their respective owners and are used only for reference.